B600

Packmaß /
Folding size

Kinnsteuerung /
Chin control

Detail Felge / Rim

Otto Bock®

QUALITY FOR LIFE

CE

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ⒺⒶ **Instructions for Use**

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Instructions for Use

Table of Contents

1	General Information.....	96			
1.1	Preface	96		4.1	Delivery
1.2	Intended Use.....	96		4.2	Initial Operation.....
1.3	Area of Application.....	97		5	Transportation and Storage.....
1.4	Service.....	98		6	Operation
2	Safety	98		6.1	Adjustment Options.....
2.1	Explanation of Symbols	98		6.1.1	Telescoping Anterior Frame
2.2	Standards and Directives	99		6.1.2	Back Rest
2.3	General Safety Information.....	99		6.1.3	Side Panel with Arm Rest.....
2.4	Safety Requirements for Transportation, Storage, and Assembly	100		6.1.4	Control Panel
2.5	Safety Requirements for Operation.....	102		6.2	Getting Into and Out of the Power Wheelchair ...
2.6	Safety Requirements for Care, Maintenance, and Disposal.....	106		6.2.1	From the Side.....
2.7	Requirements for the User	108		6.2.2	From the Front
2.8	Safety Functions.....	108		6.3	Control Unit.....
2.9	Warning Symbols and Type Plates	109		6.3.1	Control Panel
3	Product Description	112		6.3.2	Separate LCD Monitor
4	Delivery and Preparation for Use.....	113		6.3.3	Switching On and Off.....
				6.3.4	Driving Function
				6.3.5	"Battery Capacity" Indicator.....
				6.3.6	Controlling Additional Electrical Functions
				6.3.7	Drive-Away Lock

6.3.8 Lighting	132	7.7 Special Controls.....	151
6.4 Releasing and Engaging the Brake	133	7.8 Additional Indicator and Control Elements	152
6.5 Batteries	135	7.8.1 Separate LCD Monitor with Infrared.....	152
6.5.1 Charging.....	137	7.8.2 External Kilometre Indicator	154
6.5.2 Battery Charger	138	7.8.3 Push Button Module	155
7 Accessories.....	140	7.8.4 2x Adapter Cable for Piko or Buddy Button	157
7.1 Seat Options.....	141	7.9 Additional Options.....	157
7.1.1 Lift Seat Function	142	7.9.1 Lap Belt	157
7.1.2 Electric Seat Tilt	144	7.9.2 Four-Way Chest Strap.....	158
7.1.3 Electric Back Angle Adjustment.....	145	7.9.3 Adapter for Headrest Mounting Kit (Standard Seat)	159
7.1.4 Preparing Power Wheelchairs with Back Angle Adjustment for Transport	145	7.9.4 Attendant Control.....	159
7.1.5 Recaro® Seat.....	146	7.9.5 Additional Options.....	159
7.2 Caster Wheel Suspension.....	148	8 Malfunctions / Troubleshooting	160
7.3 Caster Wheel Swivel Lock	148	8.1 Warning	160
7.4 Electrically Adjustable Footrests	149	8.2 Error.....	160
7.5 Mechanically Elevating Footrest	150	8.3 Defect / Failure.....	164
7.6 Control Panel Holder	151	9 Maintenance and Care.....	164
7.6.1 Height Adjustable Control Panel Holder	151	9.1 Maintenance Intervals	164
7.6.2 Swing-Away Control Panel Holder	151	9.2 Changing the Fuse.....	167

9.3	Changing Tyres	167
9.4	Changing Defective Lights	169
9.5	Cleaning and Care	170
10	Technical Data	171
11	Disposal	173
12	Information on Re-Use	174
13	Liability	174
14	CE Conformity	175

List of Figures

Figure 1	Signage on the B600	109	Figure 20	Rear lamps	133
Figure 2	Main components.....	113	Figure 21	Brake release.....	134
Figure 3	Anti-tipper B600	114	Figure 22	Activating the brake	134
Figure 4	Fuse holder	114	Figure 23	Battery drawer	135
Figure 5	Folding size	115	Figure 24	Opening the battery drawer	136
Figure 6	Telescoping anterior frame	117	Figure 25	Snap buckle on the support strap.....	136
Figure 7	Release strap for the back rest	118	Figure 26	Batteries	136
Figure 8	Removing the side panel	118	Figure 27	Release bolt engaged	137
Figure 9	Adjusting the height of the arm rest	119	Figure 28	Charging receptacle	140
Figure 10	Adapting the control panel	119	Figure 29	Snail symbol indicating creep speed	144
Figure 11	Footrest lock	120	Figure 30	Lift seat in raised position.....	144
Figure 12	Adjusting the lower leg length.....	121	Figure 31	Electric seat tilt.....	144
Figure 13	Footplates flipped up.....	122	Figure 32	Electric back angle adjustment.....	145
Figure 14	Footrests swung to the side.....	122	Figure 33	Locking the cross bolt in the holder	145
Figure 15	Footrests removed	123	Figure 34	Illustration of pressure points.....	147
Figure 16	Control panel	124	Figure 35	Caster wheel suspension	148
Figure 17	LCD monitor with all symbols	125	Figure 36	Caster wheel swivel lock engaged and released	149
Figure 18	Information window with drive-away lock....	131	Figure 37	Elevated footrest	150
Figure 19	Front lights	132	Figure 38	Mechanically elevating footrest.....	150

Figure 39	Height adjustment of control panel holder...	151
Figure 40	Swing-away control panel holder	151
Figure 41	LCD module with infrared	152
Figure 42	Kilometre indicator	154
Figure 43	Push button module.....	156
Figure 44	Applying the lap belt.....	158
Figure 45	Adapter for headrest mounting kit	159
Figure 46	Fuse holder	167
Figure 47	Disassembling the drive wheel.....	168
Figure 48	Disassembling the caster wheel	168
Figure 49	Changing the front lights	169
Figure 50	Replacing the direction indicator lamp.....	170
Figure 51	Rear lamp, pane removed	170

1 General Information

1.1 Preface

These instructions for use provide the user as well as his or her attendants with all the required knowledge on the design, function, operation, and maintenance of the B600 power wheelchair from Otto Bock Mobility Solutions GmbH. The instructions contain all the information that is needed to use the power wheelchair safely, to determine the possible causes of a malfunction, and to help eliminate it.

Knowledge of these instructions for use is absolutely necessary for the safe operation of the power wheelchair. Therefore, the user and his or her attendants must read the instructions for use thoroughly, especially the section on "Safety", before using the power wheelchair. This will ensure that all features of the power wheelchair can be utilised to their full advantage.

1.2 Intended Use

The B600 power wheelchair is intended exclusively for individual indoor and outdoor self-transportation by persons with walking impediments or walking disabilities. The power wheelchair may only be combined with the options mentioned in these instructions for use and vice versa. Otto Bock assumes no liability for combinations with medical devices and/or accessories from other manufacturers outside of the modular system. An option for the B600 is to have it controlled by an attendant with the use of an attendant control.

Any other use is considered improper use. The manufacturer is not liable for any personal injury or damage to property resulting from improper use; in such cases, the user has sole liability.

The B600 may only be used by properly trained persons. Training in the use of the power wheelchair is required to protect persons from danger and to ensure that the B600 is operated safely and correctly.

Safe operation of the B600 can only be assured if it is used for the intended purpose in accordance with the information in these instructions for use. The user is ultimately responsible for accident-free operation.

1.3 Area of Application

The versatility of the power wheelchair and its modular design make it suitable for patients who have walking impediments or walking disabilities due to:

- Paralysis
- Loss of limb (leg amputation)
- Limb defects / deformations
- Joint contractures / defects
- Other disorders

The B600 was designed especially for users who are able to move independently in such a power wheelchair.

The following points should also be considered for custom fittings:

- Body size and body weight (max. loading capacity 140 kg, optional 180 kg)
- Physical and psychological limitations
- Age of the patient
- Living conditions
- Environment

1.4 Service

INFORMATION

Service and repairs on the power wheelchair must only be carried out by specialist dealers with personnel authorised and trained by Otto Bock. Should any problems arise, please contact your power wheelchair supplier.

Should you have any questions or a problem that cannot be resolved despite using the instructions for use, please contact Otto Bock Customer Service (see inside cover for address).

Otto Bock endeavours to support customers in all respects in order to keep them satisfied with their product for a long time.

2 Safety

2.1 Explanation of Symbols

WARNING

Warnings regarding possible risks of severe accident or injury.

CAUTION

Warnings regarding possible risks of accident or injury.

NOTICE

Warnings regarding possible technical damage.

INFORMATION

Information for using the product. Information for service personnel.

2.2 Standards and Directives

All safety information contained in these instructions for use refers to the currently valid national laws and regulations of the European Union. In other countries, compliance with the applicable laws and national regulations is required.

In addition to the safety instructions contained in these instructions for use, the user is required to observe and comply with the BGV (Employer's Liability Insurance Association regulations), UVV (accident prevention regulations), and environmental protection regulations. All information contained in these instructions for use must be complied with at all times without restrictions.

The B600 power wheelchair has been constructed in accordance with state-of-the-art technology and is safe to operate. The safety of the B600 power wheelchair is confirmed by the CE symbol and the declaration of conformity.

2.3 General Safety Information

WARNING

Risk of suffocation. Packaging materials must be kept out of reach of children.

CAUTION

Risk of accidents or injury if safety instructions are not observed or complied with. All safety instructions contained in these instructions for use and in all other applicable documents must be observed and complied with. The instructions for use must be available to the user at all times.

CAUTION

Risk of accidents or injury due to improper use. The B600 must only be used in accordance with the guidelines for proper use. The B600 must only be used by persons who have received adequate training.

CAUTION

Risk of burns in the proximity of fire. While the back rest upholstery and seat cushion of the power wheelchair are flame resistant, they may catch fire. Therefore utmost

caution near any sources of open flame or sparks, especially lit cigarettes, is required.

INFORMATION

Only use original options supplied by the manufacturer. The optional components may only be mounted as described here. Failure to comply will void the warranty.

2.4 Safety Requirements for Transportation, Storage, and Assembly

Transportation and Storage

WARNING

Risk of accidents and injury due to improper use as a seat in a vehicle for the disabled. The power wheelchair may be used in a vehicle for the disabled subject to the application of safety elements (e.g. lap belt) offered by Otto Bock and the use of suitable restraint systems. Only a single person may be transported with the B600 at one time. For optimum protection of the occupants in case of an accident, use the seats and corresponding restraint systems installed in a vehicle for the disabled while it is in operation.

INFORMATION

More information on the use of the B600 as a seat in a vehicle for the disabled is available in our brochure “Using your Wheelchair / Seating Shell Mobility Base or Buggy for Transportation in a Vehicle for the Disabled”, article no. 646D158.

NOTICE

Risk of damage due to improper transportation. Hoisting devices used for transportation must have a sufficient capacity. The power wheelchair must be secured in accordance with the requirements for the transport aid. Only attach tensioning straps to the corresponding eyelets.

During transportation on lifting platforms or in lifts, the controls of the power wheelchair must be turned off. Engage the brake.

Ensure that the power wheelchair is centred on the lifting platform.

None of the power wheelchair components, such as the anti-tipper or other parts, are allowed to be in the danger zone.

INFORMATION

Battery damage due to deep discharge. Remove the

fuse for shipping or when the power wheelchair is not being used for an extended period of time.

INFORMATION

The tyres of the power wheelchair contain chemical substances that may react with other chemical substances such as cleaning agents or acids.

Black tyres contain soot particles that can cause discolouration and black marks. A suitable mat should be used when parking the wheelchair for extended periods of time.

Assembly

⚠ CAUTION

Risk of accidents due to unsecured screw connections.

When screw connections with thread locks are loosened, they must be equipped with new thread locks or secured with a medium-strength liquid thread lock substance (e.g. EuroLock A24.20).

After the power wheelchair has been set and adjusted, the attachment screws and/or nuts must be firmly retightened. Observe the applicable torque specifications.

⚠ CAUTION

Risk of tipping due to an incorrectly mounted anti-

tipper. For safe operation, the anti-tipper must be installed correctly and maintained in proper condition.

⚠ CAUTION

Risk of accidents or injury due to improper installation of the Recaro® seat. During installation, ensure that the seat adapter is securely attached to the retainer bushings and that the anterior locking bolts have snapped into place so that they touch the key ring.

INFORMATION

Prior to using the power wheelchair, all the necessary mechanical adaptations (e. g. mounting special controls/a push button module) and software settings (e. g. programming the controls) must be made to comply with the individual requirements and abilities of the user. Settings may only be made by personnel that has been authorised and trained by Otto Bock.

INFORMATION

When installing special controls/a push button module, the user is obliged to have a specialist dealer give him or her detailed instructions on how to use it.

INFORMATION

The mounting position of the seat lift unit and / or seat must not be modified.

2.5 Safety Requirements for Operation

⚠ WARNING

Risk of accidents and injury due to defective safety functions. The user is obligated to ensure that the power wheelchair and its safety functions are in safe and proper condition before every use.

The power wheelchair may only be operated if all safety functions, e.g. the automatic brakes, are functional. Brake failure can result in serious accidents with fatal injuries.

⚠ WARNING

Risk of accidents and injury due to incorrect configuration settings. Modified parameter settings in the configuration can lead to changes in driving characteristics. In particular, changes to the speed, acceleration, braking or joystick settings can lead to unexpected and therefore uncontrollable operating performance with a risk of accidents.

Always test the driving characteristics of the power wheelchair after configuration / programming is complete. Programming must only be completed by authorised personnel. Neither Otto Bock nor the control unit manufacturer are liable for damages (especially in combination with special controls) caused by programming that was not properly / professionally adapted to the abilities of the wheelchair user.

⚠ WARNING

Risk of accidents and injury due to a lack of brake functionality when the brake is disengaged. Note the lack of brake functionality when the brake is disengaged, especially when operating the power wheelchair on slopes.

⚠ WARNING

Risk of injury due to tipping during operation. The B600 has been approved for ascending or descending inclines up to a maximum of 17 %. Navigating inclines above this percentage value is not permitted. The critical obstacle height of the B600 is 5 cm. Crossing obstacles higher than 5 cm is not permitted. Navigating stairs is not permitted.

⚠ WARNING

Risk of accidents due to tipping during operation. Reduce speed when driving downhill (select speed level 1). Obstacles must not be crossed while ascending or descending inclines. Avoid getting into or out of the wheelchair on slopes.

Before driving up slopes or over obstacles, a tilted or raised seat must be lowered to its basic position and the back be brought into the upright position. Tilting the seat back slightly is recommended when navigating downhill.

Obstacles such as steps or curbs may only be crossed at reduced speed (max. 3 km/h). Always approach obstacles at a right angle and cross over them without stopping.

⚠ WARNING

Risk of tipping when driving on unsuitable terrain. Driving on very smooth areas (e.g. icy surfaces) or extremely rough terrain (e.g. broken rock or rubble) is not permitted.

⚠ WARNING

Risk of tipping when using lifting platforms. Turn off the controls of the power wheelchair when using lifting platforms, lifts, buses, or trains. Engage the brake.

⚠ WARNING

Risk of tipping due to changes in the centre of gravity. The effects of changes in the centre of gravity on the behaviour of the power wheelchair on inclines, side hills, or when crossing obstacles must be tested with the secure support of an attendant prior to initial use.

⚠ WARNING

Risk of accidents and injury due to incorrect lifting. Attendants may only lift the wheelchair by securely mounted components, but never by the foot rests or arm rests.

⚠ WARNING

Risk of pinching in the lifting zone. When using the lift seat function, no unauthorised persons may be present in the danger zone. There must not be any interfering objects or obstacles in the lifting zone.

All attendants must be informed that there are pinch points caused by the construction in the area between the seat frame and wheelchair frame. Users and attendants must not reach into the danger zone if the electric seat tilt and / or lift seat functions are installed.

⚠ CAUTION

Risk of accident when driving without experience. Driving without experience can lead to tipping or other dangerous situations. Practice using the wheelchair on level, straightforward terrain before using it.

⚠ CAUTION

Risk of accidents due to uncontrolled driving characteristics. Uncontrolled movements can occur during the operation of the power wheelchair as a result of malfunctions. In this case, please contact your authorised dealer immediately. If any faults, defects, or other hazards that can lead to personal injury are detected, the power wheelchair must be taken out of service immediately.

⚠ CAUTION

Risk of accidents and injury due to getting into or out of the wheelchair improperly. Turn off the controls in order to get into or out of the wheelchair. The foot rests and arm rests are not capable of bearing full body weight, and therefore must not be used for getting into or out of the wheelchair.

⚠ CAUTION

Risk of accidents and injury due to rolling. Releasing the brake may result in uncontrolled rolling of the power wheelchair. Therefore, ensure that the brake is engaged after parking the power wheelchair.

⚠ CAUTION

Risk of accidents due to bad tyres. Visually inspect the tyres for sufficient tread depth and correct pressure before each use.

Avoid parking outside and in direct sunlight (UV radiation) since this causes the tyres to age quickly. As a result, the tread surface hardens and corner pieces break out of the tread. Otto Bock recommends replacing the tyres at least every 2 years regardless of the amount of wear.

If the wheelchair is parked for an extended period of time or the tyres overheat (e.g. in the vicinity of radiators or in case of exposure to strong sunlight behind glass), the tyres may become permanently deformed. Therefore, sufficient clearance should always be maintained from sources of heat, and the wheelchair should be moved frequently or placed on blocks for storage.

⚠ CAUTION

Risk of accidents due to improper clothing. Otto Bock advises users to wear light-coloured clothing or clothing with reflectors during hours of darkness.

NOTICE

Damage due to excessive heat or cold. The B600 can only be operated in the temperature range from -25 °C to +50 °C. It must not be operated at temperatures outside this range.

NOTICE

Damage caused by overload. The maximum load capacity for the B600 is 140 kg. It is possible to have the B600 retrofitted for a maximum load capacity of 180 kg.

NOTICE

Damage due to improper use of the lift seat function (optional). The lift seat function may only be used on horizontal surfaces and only with a vertical back angle. The seat function actuators are not intended for continuous use, but only for short-term limited operation (10 % use, 90 % idle time). General recommendation: Under maximum load, a pause of approximately 90 seconds is

required after 10 seconds of operation. The electric seat functions are considered independently of the driving function for this purpose.

The lift seat may only be subjected to a maximum load of 100 kg. It is possible to have the B600 retrofitted for a maximum load capacity of 180 kg.

The lift seat may not be operated in case of errors or malfunctions.

If creep speed is not activated when the lift seat is operated, please visit a specialist dealer immediately. The power wheelchair may only be used with the seat in its lowest position until the problem is corrected.

NOTICE

Interference from electromagnetic fields. The power wheelchair has been tested according to EMC regulations. The following particularities must be observed during operation:

The driving characteristics of the power wheelchair can be affected by electromagnetic fields (mobile phones or other radiating devices). Turn off all mobile devices when operating the wheelchair.

The power wheelchair can generate electromagnetic fields that can cause interference with other devices. Turn off the controls whenever you do not need them.

INFORMATION

Road traffic regulations must be observed when driving on the road.

INFORMATION

If the driving mode is changed while driving, the power wheelchair will accelerate or decelerate.

INFORMATION

The control unit of the power wheelchair is protected according to protection class IP 54 and the push button module according to protection class IP 64. Accordingly both can be used in bad weather (e. g. rain). The control unit / push button module is approved for indoor and outdoor operation and meets the requirements regarding climate and splashed water.

INFORMATION

The controls of the power wheelchair must be reactivated after an emergency stop. In the event of communication problems in the bus system of the controls, the system triggers an emergency stop and thus prevents any uncontrolled functions. If the driving function is still not available

after reactivating the controls, the brake must be released to switch over to pushing mode. In this case, contact a specialist dealer as soon as possible.

INFORMATION

If the controls do not emit an error signal as soon as the brake lever is activated, this indicates a malfunction. The settings must be checked by a specialist dealer.

INFORMATION

Manoeuvring is only allowed at reduced speeds.

2.6 Safety Requirements for Care, Maintenance, and Disposal

⚠ CAUTION

Risk of accidents and injury due to improper maintenance, repairs, or adjustments. Maintenance work on the power wheelchair may only be completed by personnel authorised and trained by Otto Bock. This also applies to all repairs and settings on the brake. Incorrect settings can lead to a loss of braking power.

⚠ CAUTION

Risk of injury due to explosive gases. Explosive gases can develop while the batteries are charging. Therefore the following safety precautions are required when charging the battery:

Turn off the controls.

Ensure sufficient ventilation when charging the batteries in an enclosed room.

Smoking and open flames are not permitted. Sparks must be avoided.

The air vents in the trim must not be covered.

⚠ CAUTION

Risk of injury due to uncontrolled movements. The fuse must always be removed for any maintenance work where the battery cover is open.

NOTICE

Unauthorized battery replacement. The battery may only be replaced by a specialist dealer. The characteristic curve of the battery charger set at the factory corresponds to the battery provided and must not be changed. Setting the characteristic curve incorrectly can result in permanent damage to the battery.

NOTICE

Risk of damage to electronics due to water penetration. Never use a water jet or high-pressure cleaning apparatus to clean the power wheelchair. Water must not come into direct contact with the electronics, motor, or battery.

NOTICE

Risk of damage due to oil leaks. If the power wheelchair is equipped with gas compression springs, watch for oil leaks on the piston rods. Leaks reduce the functional capabilities of the gas compression springs and may result in their failure. Defective gas compression springs (e. g. scratches or damage to the piston rod caused by impact; deformation of the cylinder tube) must be replaced immediately.

INFORMATION

The power wheelchair must be inspected for functional reliability and driving safety at least once a year by an authorised specialist.

INFORMATION

When closing the battery drawer, the release bolt must engage in the recess provided for this purpose. The release bolt must be easy to see and feel on the underside of the drive unit sustainer.

INFORMATION

Piston rods do not require lubrication. They are maintenance-free.

INFORMATION

Defective batteries must be disposed of properly in accordance with country-specific regulations.

2.7 Requirements for the User

⚠ CAUTION**Risk of accidents and injury due to improper operation.**

The power wheelchair may only be operated by properly informed users. For this purpose, the user and any attendants must receive instruction in the use of the power wheelchair from specialists trained and authorised by Otto Bock.

The operator must have read and understood all the information in the instructions for use.

The wheelchair may not be operated in cases of exhaustion or under the influence of alcohol or medication.

The user must not have any mental limitations which can temporarily or permanently restrict attentiveness and judgment.

2.8 Safety Functions

INFORMATION


In dangerous situations, the B600 can be turned off at any time using the on/off button. When the button is pressed, the power wheelchair brakes immediately and the electric functions will be stopped.

Malfunctions such as an insufficient supply of power to the brake are recognised by the software, triggering the emergency stop brake or reducing the speed of the power wheelchair. At the same time, a warning signal will be emitted.

2.9 Warning Symbols and Type Plates

Figure 1 Signage on the B600

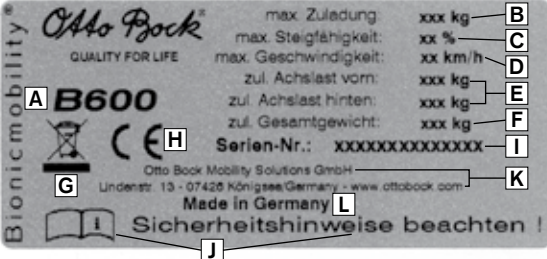Label	Meaning
	A Type designation B Maximum load capacity (see Section “Technical Data”) C Maximum climbing ability (see Section “Technical Data”) D Maximum speed (see Section “Technical Data”) E Allowable axle load front / rear F Allowable gross weight G Symbol for separate collection of electric and electronic devices. Components of the power wheelchair and batteries may not be disposed of in household waste. H CE symbol – product safety according to EU guidelines I Serial number J Read the Instructions for Use prior to using the product and observe the safety instructions mentioned there. K Manufacturer information / address L Country of origin
	Observe tyre pressure for front and rear wheels (in bar / if applicable in PSI > Pound-force per square inch)

Label	Meaning
	A Electric driving operation: lock motor brake B Manual pushing mode: unlock motor brake
	Risk of pinching. Do not reach into the danger zone.

Table 1 Signage on the B600

3 Product Description

The B600 power wheelchair is suitable for indoor and outdoor use. With a compact design, it is easy to manoeuvre indoors. The powerful drive system with two 12 V batteries combined with spring-mounted drive wheels allows obstacles to be crossed easily (Category B according to EN 12184) and offers safe driving characteristics.

The power wheelchair is controlled by an enAble50 wheelchair control system. It includes a control panel to enter driving commands and display the current status as well as a controller that controls the drive motors and other electrical functions based on the input data. Data transmission is realised via a bus system.

The enAble50 is programmable, which means the control system can be customised to individual user requirements. For example the speed, acceleration, and delay settings can be adjusted.

Special features of the B600 include:

- Adjustable centre-of-gravity shifting by means of a telescoping frame.
- Serviceability due to easy, straightforward access to all components.
- Customisation with options and custom fabrication using modular components.

The modular design of the B600 makes it possible to add other modules and devices in addition to the main components (see Figure 2), e. g. electric seat adjustment, special controls, tray, etc. These options are described in Section 7.

4 Delivery and Preparation for Use

4.1 Delivery

INFORMATION

The options included in the scope of delivery depend on the product configuration that has been selected for the power wheelchair.

The scope of delivery includes:

- Adapted power wheelchair with main components (see Figure 2)
- Battery charger
- Instructions for use
- Options (see Section 7)

The specialist dealer delivers the power wheelchair ready for use. All settings correspond to the indications on the order form or are adjusted directly on site by the dealer. The power wheelchair is adapted to the personal requirements of the individual.

The functionality of the individual components can be tested by following the instructions in Section 6. Section 8 includes a description of possible malfunctions.

Figure 2 Main components

- | | |
|-----------------------------|--------------------------|
| 1 Back rest | 5 Footrest |
| 2 Joystick and control unit | 6 Motor with drive wheel |
| 3 Arm rest (side panel) | 7 Anti-tipper |
| 4 Seat cushion | |

Figure 3 Anti-tipper B600

4.2 Initial Operation

⚠ WARNING

Risk of suffocation. Packaging materials must be kept out of reach of children.

All components must be checked for completeness and functionality (see Figure 2) prior to initial operation. Before the power wheelchair can be switched on, the 80 A fuse must be inserted into the fuse holder at the rear of the battery case (see Figure 4).

To do so, remove the fuse from the protective cover on the control panel, open the cover of the fuse holder, and insert the fuse.

Make sure that the fuse is pressed into the centre of the spring contacts and that it is not positioned at an angle. Close the cap until you feel it snap back into place.

Figure 4 Fuse holder

- 1 Cover open
- 2 Fuse inserted
- 3 Fuse holder

5 Transportation and Storage

⚠ WARNING

Risk of injury due to improper transportation. Pre-installed vehicle seats with corresponding restraint systems offer optimum protection for occupants in case of an accident. The B600 power wheelchair should only be used as a seat in a vehicle for the disabled in exceptional cases, and only in conjunction with the safety elements and restraint systems offered by Otto Bock. For more information, please refer to our brochure "Using your Wheelchair / Seating Shell Mobility Base or Buggy for Transportation in a Vehicle for the Disabled", article no. 646D158.

⚠ CAUTION

Risk of injury due to insufficient restraints. Secure the power wheelchair during transportation in another vehicle by tightening the tensioning straps properly.

The frame of the power wheelchair has 4 eyelets for the attachment of tensioning straps.

Before transporting the power wheelchair, switch off the controls and engage the brake.

You can reduce the size (Figure 5) of the power wheelchair for transportation by folding down the back rest and removing the side panels and footrests (see Section 6.1).

Figure 5 Folding size

Tyres contain chemical substances that can react to other chemical substances (such as cleaning agents, acids, etc.).

Black tyres also contain soot particles.

The power wheelchair must be stored in a dry place.

For transportation and storage, an ambient temperature between -10 °C and +40 °C must be observed.

INFORMATION If your power wheelchair is not moved for several days, permanent colour changes may occur where the wheelchair comes into contact with the surface it is standing on. A suitable mat should be used when parking the wheelchair for extended periods of time.

INFORMATION Black tyres may leave black marks on the ground under some circumstances. Therefore we recommend choosing grey tyres if the power wheelchair is mainly used indoors.

INFORMATION Direct exposure to sunlight / UV radiation causes the tyres to age more quickly. As a result, the tread surface hardens and corner pieces break out of the tread.

INFORMATION Avoid parking the wheelchair outdoors whenever possible. The tyres should be replaced at least every 2 years regardless of wear and tear.

INFORMATION Remove the fuse for shipping or when the power wheelchair is not being used for an extended period of time.

6 Operation

6.1 Adjustment Options

CAUTION

Risk of accidents due to unsecured screw connections

When screw connections with thread locks are loosened, they must be equipped with new thread locks or secured with a medium-strength liquid thread lock substance (e.g. EuroLock A24.20).

After the power wheelchair has been set and adjusted, the attachment screws and/or nuts must be firmly retightened. Observe the applicable torque specifications.

Various adjustments can be made to the power wheelchair.

The seat height, seat width, and seat angle have been set in accordance with the indications on the customer order form and may only be changed by the specialist dealer.

The rear wheel suspension, which is included in the standard equipment, must be adjusted by an authorised dealer.

The following items can be adapted by the user:

- Back angle
- Arm rest height
- Arm rest position
- Lower leg length

The footrests and side panels can be removed if needed.

6.1.1 Telescoping Anterior Frame

The telescoping anterior frame of the power wheelchair allows the wheelbase to be adapted to the respective seat (Figure 6). Since this setting involves a change to the centre of gravity, it may only be carried out by an authorised dealer.

Figure 6 Telescoping anterior frame

6.1.2 Back Rest

The back rest can be set to 4 different positions by pulling the release strap (base of the back rest, Figure 7).

Once the back rest is in the desired position, release the strap to engage the locks.

Figure 7 Release strap for the back rest

6.1.3 Side Panel with Arm Rest

The side panels can be removed by loosening the thumb screws and lifting the side panels. If the control panel is attached to the arm rest, it must be removed before taking out the side panel (Section 6.1.4).

The thumb screws are located at the lower end of the arm rest holders (see Figure 9).

Figure 8 Removing the side panel

After remounting the side panel, do not forget to firmly retighten the thumb screw.

The height of the arm rest can be changed by loosening the screw (see Figure 9).

Firmly re-tighten the screws after making adjustments.

Figure 9 Adjusting the height of the arm rest

6.1.4 Control Panel

NOTICE

Cable damage. Positioning the cable incorrectly can lead to pinching and damage to the cable. The cable must not be too tight or too loose. Avoid bending or squeezing the cable.

Adapting the Control Panel to Arm Length

To adapt the control panel to the arm length, the 3 screws on the bottom of the arm rest must be loosened with an Allen key (size 3) (see Figure 10). Now the control panel can be moved backwards and forwards. The arm rest can be adjusted in the same way.

Figure 10 Adapting the control panel

Firmly re-tighten the screws after making adjustments.

INFORMATION If the control panel rail is too long, the protruding part can simply be cut off with a saw.

Mounting and Removing the Control Panel

The standard configuration has the control panel mounted on the right side. If required, it is also possible to mount the control panel on the left side. This change may only be made by a specialist dealer.

6.1.5 Footrest

⚠ CAUTION

Risk of pinching. Ensure fingers are not in the danger area when flipping the footplates up or down. Complete removal of the footrest from the adapter is only permitted for transportation of the wheelchair.

Removal

1. Remove the calf band.
2. Flip the footplate up, release the footrest lock (see Figure 11), and swing the footrest either to the inside or outside.
3. Pull the footrest up to remove it.

Reassembly

1. Reinsert the footrest from above by placing it in the adapter and swinging it forward until the footrest lock engages.
2. Reattach the calf band to the adapter.

Figure 11 Footrest lock

Adjusting the Lower Leg Length (see Figure 12)

By loosening the set screw, the footplate can be adapted to your lower leg length and to the thickness of your seat cushion.

INFORMATION When making the adjustment, ensure that the footrest bar (item b) slides at least 60 mm into the swivel segment (item c). After completing the adjustment, tighten the set screw with a torque of 8 Nm.

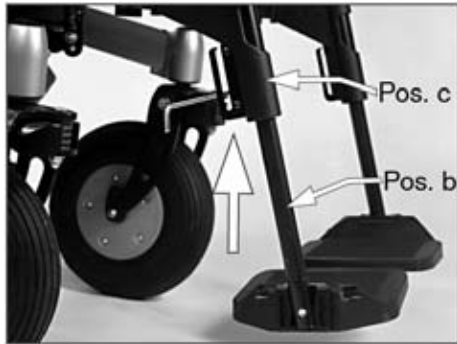
Figure 12 Adjusting the lower leg length

6.2 Getting Into and Out of the Power Wheelchair

⚠ CAUTION

Risk of injury if the power wheelchair starts rolling. When getting into or out of the power wheelchair, the wheelchair controls must always be switched off. This will automatically engage the motor brake.

⚠ CAUTION

Risk of breakage due to overloading. The user must not apply his or her entire body weight to the footrests and arm rests for support when getting into or out of the wheelchair.

The modular design of the power wheelchair and the ease with which you can remove the side panels and footrests make it easy to get into and out of the wheelchair from the side or from the front.

Users can get into and out of the wheelchair using their preferred and individual method.

6.2.1 From the Side

Getting into the wheelchair from the side requires the user to remove the right or left side panel depending on the side

where the user will get into the wheelchair. Bring the power wheelchair as close as possible to the place where its user is sitting.

If the control panel is located on the side where the user is getting in or out, carefully remove the control panel. Remove the side panel as described in Section 6.1.3 (remove the footrest if necessary; see Section 6.1.5).

The user can then slide onto the wheelchair seat from the side. Use of a transfer board will make this even easier.

6.2.2 From the Front

Flipping up the two footplates makes it possible for the user to get into and out of the wheelchair from the front (Figure 13). Swinging the footrests away to the sides increases the amount of space available when getting in or out (Figure 14). If required, it is also possible to remove the footrests from their adapter after swinging them out of the way (Figure 15; please observe the safety information in Section 6.1.5, footrest).

The assistance of an attendant or a transfer lifter allows the user to get into or out of the power wheelchair easily. Use of a rotation plate is also possible.

Figure 13 Footplates flipped up

Figure 14 Footrests swung to the side

Figure 15 Footrests removed

6.3 Control Unit

NOTICE

Risk of impairing the driving performance of the power wheelchair. The driving characteristics of the power wheelchair can be affected by electromagnetic fields (mobile phones or other radiating devices). Therefore all mobile devices must be turned off when driving.

NOTICE

Risk of damage to other devices. The power wheelchair can generate electromagnetic fields that can cause interference with other devices. Therefore, switch off the controls whenever you do not need them.

6.3.1 Control Panel

The power wheelchair is operated using the control panel.

The control panel consists of a keypad, LCD monitor, and joystick. The programming receptacle and two inputs for external push-buttons are located on the underside. The power wheelchair is turned on and off using the keypad; driving commands can be entered and the status of certain functions and components can be displayed.

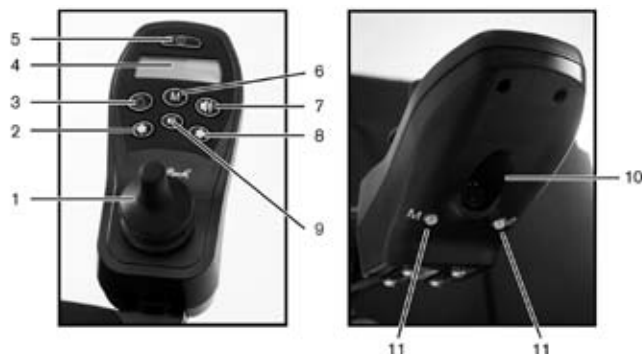
Figure 16 Control panel

- | | |
|-----------------------------|---------------------------------------|
| 1 Joystick | 6 Mode button |
| 2 Direction indicator, left | 7 Horn |
| 3 Warning flasher | 8 Direction indicator, right |
| 4 LCD monitor | 9 Light |
| 5 On/Off button | 10 Programming receptacle (underside) |
| | 11 Inputs for external push-buttons |

On/Off Button

The on/off button is used to turn the power wheelchair on, activate the drive-away lock, and turn the power wheelchair off.

Mode Button

Briefly pressing the mode button increases the speed level. After reaching the maximum speed level, you can change back to speed level 1 by pressing the button again. Pressing and holding the button (at least 2 seconds) activates the “Additional electrical functions” menu (see Section 6.3.6).

Horn

The horn will sound as long as the horn button is being pressed.

Warning Flasher

Pressing the warning flasher button activates / deactivates all four direction indicator lights.

Light

The front and rear lights are activated and deactivated by pressing the light button.

Direction Indicator Lights

The right or left front and rear direction indicator lights are activated and deactivated by pressing the “direction indicator light, right” or “direction indicator light, left” buttons respectively. The direction indicator lights automatically turn off after 20 seconds.

LCD Monitor

The LCD monitor is the communication interface between the user and the control unit. It indicates the selected speed level, the remaining battery capacity, and the status of electric options and special functions as well as warnings and errors. All display symbols are visible during the initialisation phase.

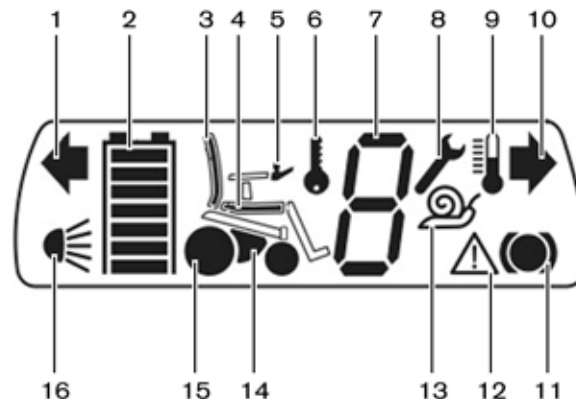
Figure 17 LCD monitor with all symbols

- | | |
|------------------------------------|--------------------------------------|
| 1 Direction indicator lights, left | 9 Excess temperature |
| 2 Battery capacity | 10 Direction indicator lights, right |
| 3 Electric back rest | 11 Brake for drive wheels |
| 4 Electric seat tilt | 12 Warning |
| 5 Control panel | 13 Creep speed |
| 6 Drive-away lock | 14 Power module |
| 7 Speed level | 15 Drive motor |
| 8 Open end wrench | 16 Light |

6.3.2 Separate LCD Monitor

A separate LCD monitor with infrared and Bluetooth (optional) is available to operate special controls and for environment control functions. Please see Section 7.8.1 for more information.

6.3.3 Switching On and Off

WARNING

Risk of death due to lack of brake functionality. The brake release lever must be locked when using the power wheelchair. The automatic brake must be operational and functional.

CAUTION

Risk of accidents due to incorrect tyre pressure. The power wheelchair must be visually inspected prior to each use to ensure that the tread depth is sufficient and that the tyres are inflated to the correct pressure. Incorrect tyre pressure reduces the life of the tyres and impairs the driving characteristics.

The control unit of the power wheelchair is switched on and off by pressing the on/off button (see Figure 16, item 5).

If the control is not used for a certain time, which is adjustable, the wheelchair is switched off automatically. It is also possible to switch off the power wheelchair with the on/off button while driving. In this case, the wheelchair brakes immediately until it stops.

INFORMATION Every time you switch on the control unit, it will return to the previously selected speed level (factory setting; can be adapted by the specialist dealer).

6.3.4 Driving Function

CAUTION

Risk of injury due to uncontrolled operation of the power wheelchair. Switch off the control unit on the power wheelchair whenever you do not need it. This will prevent you from unintentionally activating the joystick.

CAUTION

Risk of injury due to tipping during operation. Observe the following safety information when operating the power wheelchair:

- Operation of the wheelchair is limited to inclines up to a maximum of 17 %.

- When driving downhill, reduce speed in accordance with the degree of the incline.
- Do not cross obstacles higher than 5 cm.
- Do not cross steps or curbs without reducing speed.

The power wheelchair has been approved for ascending or descending inclines of up to a maximum of 17 %. Driving on steeper inclines or slopes is not permitted.

In order to navigate downhill grades safely, the speed must be reduced according to the slope (e.g. select speed level 1).

The critical obstacle height of the power wheelchair is 5 cm. Crossing obstacles higher than 5 cm is not permitted. Obstacles such as steps or curbs must be crossed at reduced speeds.

If there are obstacles in the travel path, it is important to drive around them with a large margin of safety. On uneven ground, the driving characteristics of the wheelchair may get out of control. Therefore the speed must always be adapted to the ground conditions.

The battery indicator and speed level are shown in the driving menu on the LCD monitor.

The joystick (see Figure 16, item 1) is used for driving. The further the joystick is moved away from the centre position, the faster the power wheelchair will drive in this direction. The maximum speed with full deflection of the joystick depends on the selected speed level. Releasing the joystick automatically activates the brake function, which brings the wheelchair to a halt. When standing still, the mechanical brake are automatically active so the power wheelchair cannot roll.

The power wheelchair has a programmable number of speed levels (delivery condition = 5 speed levels). The "M" button (see Figure 16, Item 6) is used to increase the speed level. After reaching the maximum speed level, the control unit switches back to speed level 1.

Adapting the Driving Characteristics

WARNING

Risk of accidents and injury due to incorrect configuration settings. Modified parameter settings in the configuration can lead to changes in driving characteristics. In particular, changes to the speed, acceleration, braking or joystick settings can lead to unexpected and therefore uncontrollable operating performance with a risk of accidents.

Always test the driving characteristics of the power wheelchair after configuration / programming is complete.

Programming must only be completed by authorised personnel. Neither Otto Bock nor the control unit manufacturer are liable for damages (especially in combination with special controls) caused by programming that was not properly / professionally adapted to the abilities of the wheelchair user.

The speed, acceleration, and delay values can only be adjusted to individual user requirements by the specialist dealer via a hand-held programming device that is connected to the control panel.

6.3.5 “Battery Capacity” Indicator

INFORMATION

Immediately after switching the power wheelchair on, the battery indicator shows the battery capacity that was saved before the wheelchair was switched off the last time. The exact battery status is displayed after driving a short distance.

The battery indicator on the LCD monitor is divided into 7 segments that show the current battery capacity.

A charge of 100% corresponds to 7 segments on the battery symbol. Every time one segment disappears, this means that the battery capacity has decreased by approximately 14%. At a total range of approximately 35 km over even terrain, each segment corresponds to a range of approximately 5 km when the speed is kept constant.

When the last segment flashes, the battery must be charged immediately.

If all the segments have disappeared and the battery symbol is flashing, the battery is under voltage. Since any further use will lead to battery damage, a warning signal also appears.

If the battery symbol and all the segments are flashing, this means that the battery is over voltage. Since this will lead to battery damage, a warning symbol also appears.

The charging process is indicated by the battery segments lighting up one after the other. When the battery is charging, the driving function is locked out (see Section 6.5.1, Charging).

Display symbol(s)	Information
	Driving menu with speed level and full battery capacity
	Low battery capacity
	Charging process with drive function lock
	Battery is under voltage, with warning symbol
	Battery is over voltage, with warning symbol

Table. 2 Battery indicator on the control panel

6.3.6 Controlling Additional Electrical Functions

INFORMATION

For more information on additional electrical functions, please see Section 8.1 through 8.3.

Additional electrical functions are accessed by pressing the “M” button (approximately 2 seconds). Moving the joystick to the right lets you switch between the different functions. The function that is currently selected is shown on the LCD monitor.

The respective function can be extended or retracted by moving the joystick ahead or back. The electric motor moves the respective function as long as the joystick is deflected, and stops at the end positions.

Electric back angle adjustment: When the joystick is deflected ahead or back, the back rest moves ahead or back accordingly.

Electric seat tilt: Deflecting the joystick to the rear will electrically tilt the seat to the rear. Deflecting the joystick to the front tilts the seat ahead until it is level (factory setting; can be adapted by the specialist dealer).

Lift seat: Deflecting the joystick ahead or back raises or lowers the electric lift seat. As soon as the lift seat leaves the lowest position, the speed is reduced (creep speed) because the power wheelchair is less tip-resistant. This is indicated by a snail symbol on the control panel LCD monitor. (see Figure 17).

Electric footrests: Deflecting the joystick ahead or back raises or lowers the electric footrests.

Coupled seat adjustment: Deflecting the joystick ahead or back moves the back rest and seat forward and down or back and up simultaneously.

You can switch back to driving mode by pressing the “M” button.

Display symbol(s)	Information
	Electric back angle adjustment
	Electric seat tilt
	Lift seat
	Electric footrests
	Coupled electric back angle adjustment and seat tilt

Table 3 Display of electric seat functions on the control panel

6.3.7 Drive-Away Lock

INFORMATION

The drive-away lock can be activated as described in the following, only if the “drive-away lock” parameter has been set to “On”. At the factory, the drive-away lock parameter is set to “Off”. It is also possible that the factory setting has been changed to “On” if requested with the order or by the specialist dealer.

Please ask your specialist dealer about the setting selected for your wheelchair.

The control unit of the B600 power wheelchair features an electric drive-away lock to disable the driving function.

Activation via control panel

1. If the control unit is turned on, press and hold the mode button for at least 5 seconds.
2. A short beep confirms activation of the drive-away lock.
3. The control unit is turned off automatically.

Activation of the drive-away lock is indicated by the key symbol on the display.

Display symbol(s)	Information
	Drive-away lock

Table 4 Display of the drive-away lock on the control panel

Activation via separate LCD monitor

If a special control is used, the drive-away lock is activated using the separate LCD monitor. On the LDC monitor, select the menu item “Drive-away lock” in the menu “Settings”. Activation is confirmed with an acoustic signal and the control unit is switched off.

Figure 18 Information window with drive-away lock

Deactivation via control panel

1. Press the On/Off button to turn on the control unit.
2. Move the joystick to the front end position until you hear a beep.
3. Move the joystick to the rear end position until you hear a beep.
4. Release the joystick.
5. The speed level and battery indicator appear on the LCD monitor, and the control unit is ready for driving.

Now the drive-away lock is deactivated and the wheelchair can be driven.

INFORMATION If the joystick is not moved correctly, the drive-away lock remains activated. In order to deactivate the drive-away lock, the control unit must be switched off. The power wheelchair can be switched on again and the drive-away lock can be deactivated.

Deactivation via separate LCD monitor

After switching the power wheelchair on, an information window appears on the LCD monitor (see Figure 18). The deactivation process is the same as described under “Deactivation via control panel”.

6.3.8 Lighting

The B600 power wheelchair is equipped with lighting as standard equipment. One halogen light and one direction indicator is attached to each side panel on the right and left (see Figure 19).

Figure 19 Front lights

Two back-up lamps with integrated direction indicator lights are integrated into the rear covering (see Figure 20). The control panel is used to turn on the warning flashers, lights, and right/left direction indicator lights (see Figure 16).

Figure 20 Rear lamps

6.4 Releasing and Engaging the Brake

⚠ WARNING

Risk of death due to brake malfunction. Incorrect brake settings can lead to a loss of brake functionality and therefore to serious bodily injuries and even death. Repairs and adjustments to the brake may only be carried out by authorised service personnel.

⚠ WARNING

Risk of accidents and injury due to lack of brake functionality. When the brake is deactivated (push mode), no brake functionality is available. When moving the power

wheelchair on an incline, the person pushing must provide the appropriate brake force.

NOTICE

Risk of damage when parking the wheelchair with the brake deactivated. Deactivating the brake may result in uncontrolled rolling of the power wheelchair. Therefore, make sure that the brake is engaged when parking the power wheelchair.

In case of control unit failure or insufficient battery capacity, the power wheelchair can be pushed. To do so, the brake is released by mechanically unlocking it. The brake release mechanism is located between the frame and seat.

Releasing the Brake

For safety reasons, you first need to pull up on the brake release bolt (see Figure 21).

INFORMATION In this position, the brake still has to be engaged! If the pushing function is already active, please contact your specialist dealer.

INFORMATION When the joystick is activated while the brake is deactivated, the control unit emits an error signal on the control panel. If this is not the case, there is a malfunction that must be corrected immediately by a specialist dealer.

**Figure 21 Brake release**

- 1 Brake release bolt 2 Brake release lever

Now the brake release lever has been released and can be pushed ahead to the stop. In this position, the control unit will recognise that the brake has been released and automatically deactivate the driving function. As soon as the joystick is deflected, a warning will appear on the LCD monitor.

INFORMATION Once the brake release lever has been released, all braking systems are deactivated.

Activating the brake

To activate the brake, push the brake release lever up until the brake release bolt engages (see Figure 22). To reactivate the driving function, switch the wheelchair control unit off and then back on.

**Figure 22 Activating the brake**

6.5 Batteries

INFORMATION

When checking the power wheelchair's batteries, make sure it is standing on a level surface.

INFORMATION

Before doing any maintenance work on the batteries, please read the battery manufacturer's warnings thoroughly!

INFORMATION

The fluid level in the batteries has to be checked monthly. Top up the batteries with distilled water when required.

The standard version of the B600 power wheelchair includes two low-maintenance 12 V wet cell batteries.

The batteries are located under the seat of the power wheelchair.

For maintenance work and for mounting/removing the batteries, the B600 power wheelchair has a battery drawer (Figure 23).

Figure 23 Battery drawer

To check the batteries, proceed as follows:

1. Push the release bolt of the battery drawer up. The release bolt is located in a recess in the rear covering (see Figure 24).
2. Pull out the battery drawer towards the rear.
3. Undo the snap buckle on the support strap for the battery drawer (see Figure 25) and remove the cover. Once the cover has been removed, you have unrestricted access to the batteries (see Figure 26).

4. Unscrew the caps (see Figure 26) of the individual chambers with a large screwdriver. You will be able to see the battery fluid level on the acid level markings.
5. If the fluid is below this marking, top up with distilled water.
6. Close the battery case.
7. Close the battery drawer.

Maintenance-free lead gel batteries can be used as an alternative.

Figure 24 Opening the battery drawer

Figure 25 Snap buckle on the support strap

Figure 26 Batteries

INFORMATION Ensure the release bolt engages properly when closing the battery drawer. It has to lock into the corresponding recess on the drive unit sustainer.

When it is locked correctly, the bolt is easy to see and feel on the underside of the drive unit sustainer (see Figure 27).

Figure 27 Release bolt engaged

6.5.1 Charging

⚠ WARNING

Risk of injury due to explosive gases. Explosive gases can develop while the batteries are charging. The following safety instructions must be followed under all circumstances: Ensure sufficient ventilation when charging the batteries in an enclosed room. Smoking and fire are not permitted. Sparks must be avoided. Do not cover the air vents in the trim.

NOTICE

Risk of battery damage. Driving with low battery capacity for an extended period of time will discharge the battery completely and result in damage. There is a risk that the power wheelchair may stop due to zero battery capacity and bring the user into a dangerous situation.

The remaining battery capacity determines the range of the power wheelchair. The battery capacity is influenced by many factors. In addition to the temperature, age of the batteries, and amount of use, the charging cycle has a pronounced effect on the capacity and therefore on the range.

For an optimum charging frequency, the following applies:

- The batteries can be charged at any time, regardless of the remaining charge.
- It takes about 10 hours until a discharged battery (only one flashing segment) is fully charged. When the charging process is complete, the battery charger can remain connected to the power wheelchair with no risk of overcharging or damaging the battery. The battery charger features a programmed recharging phase that will maintain the battery capacity at the level that has been reached.
- For daily use, connecting the charger overnight is recommended so that the full capacity is available every day.
- If the wheelchair is not used for extended periods of time, the batteries gradually discharge. If the power wheelchair is not driven for extended periods of time, a charging cycle must be carried out at least 1 x per week in order to maintain the battery capacity. Removing the fuse is recommended for extended periods of disuse.
- Total discharge of the batteries should be avoided.
- When the batteries are charging, the power wheelchair control unit must be switched off to allow all of the charging current to be fed into the battery.

The following must be observed when charging the batteries:

- Only the battery charger provided by Otto Bock for this purpose may be used for charging. Failure to comply with this requirement will render the warranty null and void.
- The voltage settings on the battery charger must correspond to the voltage used in your country.

6.5.2 Battery Charger

WARNING

Explosion hazard due to sparks. Always switch off the battery charger and remove the mains plug before disconnecting the battery.

NOTICE

Unauthorized battery replacement. The battery may only be replaced by a specialist dealer. The characteristic curve of the battery charger set at the factory corresponds to the battery provided and must not be changed. Setting the characteristic curve incorrectly can result in permanent damage to the battery.

The battery charger is designed for maintenance-free and low-maintenance batteries. The charger can operate along two characteristic curves to charge the respective battery types. The power wheelchair is delivered with the appropriate characteristic curve set in accordance with the battery type. If the battery charger is to be used with another power wheelchair or if new batteries are installed, the setting has to be checked.

The battery type can be selected by means of a switch located on the bottom of the battery charger (covered by a green adhesive label). It is indicated by the green flashing LED at the beginning of the charging process.

NOTICE **Risk of battery damage.** Setting the characteristic curve incorrectly can result in permanent damage to the battery.

When handling the battery charger, the following safety instructions must be observed. Failure to do so can impair the functionality of the battery charger:

- Always place the rubber feet of the battery charger on a level surface.
- Protect the battery charger from direct sunlight to avoid additional heat.

- The location where the battery charger is set up must be dry and well ventilated. Avoid the penetration of dust and dirt.

- Use a dry cloth to clean the battery charger.

To charge the batteries, proceed as follows:

1. Turn off the control unit on the power wheelchair.
2. Open the hinged flap of the charging receptacle; the charging receptacle is integrated in the left or right lateral covering (see Figure 28).
3. Insert the battery charger plug into the charging receptacle on the power wheelchair.
4. Connect the battery charger to a wall outlet and turn it on. The batteries will start charging automatically, and the current charging status is indicated on the LEDs of the battery charger (see Table 5).
5. When the charging process is complete, turn off the battery charger: Disconnect the battery charger from the mains supply and then from the power wheelchair.
6. Close the charging receptacle cover to protect the receptacle against moisture.
7. Turn on the control unit; the power wheelchair is now ready for operation.

Figure 28 Charging receptacle

The following statuses are indicated by the battery charger:

Display	Function
Green LED flashes 1x	Low-maintenance battery
2x	Maintenance-free battery (gel)
steady	Battery is charged to 80 %
Yellow LED is illuminated	Battery is charging
Green LED is illuminated	Battery is fully charged
Red LED is illuminated	Battery is defective or missing, or charging time has been exceeded

Table 5 Battery charging status during the charging process

If none of the LEDs are lit, there is no power supply.

7 Accessories

⚠ CAUTION

Risk of accidents due to unsecured screw connections.

When screw connections with thread locks are loosened, they must be equipped with new thread locks or secured with a medium-strength liquid thread lock substance (e.g. EuroLock A24.20).

After the power wheelchair has been set and adjusted, the attachment screws and/or nuts must be firmly retightened. Observe the applicable torque specifications.

INFORMATION

Only original options provided by the manufacturer may be used. The optional components may only be mounted as described here. Failure to comply will void the warranty.

INFORMATION

All available optional add-on components are contained on the order form and in the wheelchair accessories catalogue.

The B600 power wheelchair has been designed as a modular system. Certain component groups can be exchanged and other accessories can be added.

The enAble50 control unit facilitates the use of numerous additional electrical functions as well as customised input devices.

The complete range of options is listed on the order form and in the wheelchair accessories catalogue.

7.1 Seat Options

WARNING

Risk of injury due to tipping during operation. Before driving up slopes or over obstacles, a tilted seat must be lowered to its basic position, the lift seat must be in its lowest position, and the back must be in the upright position. Slightly tilt the seat to the rear when driving down inclines or over obstacles and reduce the speed (max. 3 km/h).

WARNING

Risk of pinching or crushing fingers. Do not reach into the danger zone between the frame and the seat when

tilting the seat. There must not be any interfering objects or obstacles in the tilting zone.

CAUTION

Risk of injury when tilting the back rest. Overloading the actuators can result in spindle nut breakage and thus cause the seat to drop or the back rest to fold back.

INFORMATION

It is easier to tilt the seat if the back rest is inclined to the rear.

INFORMATION

The seat function actuators are not intended for continuous use, but only for short-term limited operation. General recommendation: At maximum load capacity, 10 seconds of operation should be followed by 90 seconds of idle time. The electric seat functions are considered independently of the driving function for this purpose.

INFORMATION

When getting into or out of the power wheelchair, bring the seat to a level position.

7.1.1 Lift Seat Function

CAUTION

Risk of injury if the power wheelchair tips over. The power wheelchair can tip over when the lift seat is raised, causing the user to fall out of the wheelchair. The following safety instructions must be followed under all circumstances:

- Only travel short distances when the lift seat is raised, e. g. when manoeuvring.
- Use the lift seat function only with the back rest set to a vertical position.
- Apply the lap belt if possible.
- Only use the lift seat function on level surfaces!
- Do not change the mounting position of the lift seat mechanism or the lift seat.
- For lift seat retrofit or conversion, observe the instructions provided in the service manual.

WARNING

Risk of injury if the power wheelchair tips over. Creep speed is automatically activated when the lift seat is operated. If this is not case, it represents a malfunction and the

lift seat function must not be used. Contact a specialist dealer as soon as possible to eliminate the problem.

WARNING

Risk of pinching in the lifting zone. When using the lift seat function, no unauthorised persons may be present in the danger zone. There must not be any interfering objects or obstacles in the lifting zone.

All attendants must be informed that there are pinch points caused by the construction in the area between the seat frame and wheelchair frame. Users and attendants must not reach into the danger zone if the electric seat tilt and / or lift seat functions are installed.

CAUTION

Risk of falling. Do not lean the upper body beyond the edge of the seat when the lift seat is raised.

CAUTION

Risk of accidents due to a restricted field of vision. Since the user's field of vision is restricted when the seat is raised, extra caution is required when operating the wheelchair with the lift seat in the raised position and when activating the lift seat function.

⚠ CAUTION

Risk of tipping in case of additional use of the seat tilt function. The seat tilt function must be set to the horizontal position when driving up inclines or over curbs. Since tipping resistance is reduced in such situations, additional caution is required!

NOTICE

Risk of damage due to improper maintenance. The lift seat must be checked for visible damage and the tightness of all screw connections must be verified at least once a month. Also maintain correct air pressure.

NOTICE

Risk of damage during transportation. Always lower the lift seat fully for loading or transportation.

NOTICE

Risk of damage due to overloading. The maximum load capacity of the lift seat is 100 kg for the standard version or 180 kg for the heavy-duty version.

INFORMATION

A sufficient supply of power must be maintained in order to guarantee proper operation of the lift seat.

INFORMATION

The lift seat is not intended for continuous use, but only for short-term limited operation (10 % use, 90 % idle time). For every 1 hour of operation, the lift seat function can be used for a maximum of 6 minutes.

The B600 can be equipped with an optional lift seat. With the lift seat function, the seat can be raised by up to 40 cm using a motor (Figure 30). It is controlled by the joystick in “additional function” mode or by the push button module (see Figure 41).

The drive function remains available even when the lift seat is raised. As soon as the lift seat leaves the lowest position, the speed is reduced (creep speed) because the power wheelchair is less tip-resistant. This is indicated by a snail symbol on the control panel or the optional, separate LCD monitor (Figure 29).

Figure 29 Snail symbol indicating creep speed

Figure 30 Lift seat in raised position

7.1.2 Electric Seat Tilt

INFORMATION

For safe operation, please also observe the information in Section 2.3 “General Safety Information” and Section 2.5 “Safety Requirements for Operation”.

The electric seat tilt adjustment lets you bring the seat into a tilted position of up to 30°, e.g. for pressure relief. It is controlled by the joystick in “additional function” mode or by the push button module. The seat can be tilted backwards continuously (Figure 31). The power wheelchair has an integrated feature that shifts the centre of gravity and serves to improve stability against tipping.

Figure 31 Electric seat tilt

7.1.3 Electric Back Angle Adjustment

INFORMATION

For safe operation, please also observe the information in Section 2.3 "General Safety Information" and Section 2.5 "Safety Requirements for Operation".

The back rest can be equipped with an electric back angle adjustment option (see Figure 32). It is controlled by the joystick in "additional function" mode or by the push button module.

The backrest can be inclined backwards continuously by up to 30°.

Figure 32 Electric back angle adjustment

7.1.4 Preparing Power Wheelchairs with Back Angle Adjustment for Transport

To obtain the folding size, the following must be done with the back angle adjustment:

- Release the cross bolt at the lower end of the drive unit (by pushing up the lever, see arrow in Figure 33).
- Remove the bolt from the holder.
- Remove the side panels and fold down the seat back onto the seat bottom.

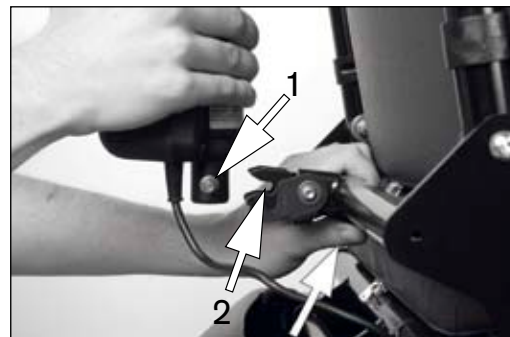
Figure 33 Locking the cross bolt in the holder

- | | |
|--------------|-------------------|
| 1 Cross bolt | 2 Lock with lever |
|--------------|-------------------|

To put a power wheelchair with back angle adjustment into operation:

- Unfold the back rest.
- Insert the cross bolt at the end of the drive unit into the holder.
- Lock the bolt (by pushing down the lever) so the latch engages the bolt (Figure 33).

INFORMATION Verify that the back angle adjustment bolt and lock are secure!

7.1.5 Recaro® Seat

⚠ WARNING

Risk of injury due to tipping during operation. Before driving up slopes or over obstacles, a tilted seat must be lowered to its basic position, the lift seat must be in its lowest position, and the back must be in the upright position. Slightly tilt the seat to the rear when driving down inclines or over obstacles and reduce the speed (max. 3 km/h).

⚠ CAUTION

Risk of tipping due to incorrect back angle adjustment. When standing still, the back rest may be inclined to the rear up to 30° only. During driving on level ground, the back rest may be inclined to the rear by up to 20° only.

INFORMATION

For safe operation, please also observe the information in Section 2.3 “General Safety Information” and Section 2.5 “Safety Requirements for Operation”.

The power wheelchair can be equipped with various Recaro® seat models.

To adjust the back angle of a Recaro® seat model, use the rotary handle on the right or left side of the back rest.

A release handle is used to fold down the back rest. It is located on the side of the back rest. By pulling the release handle upward, the back rest can be folded down to the front. The back rest automatically locks in place when it is brought into its initial position. Verify that the back rest is securely locked in place every time it has been unlocked.

The Recaro® seat can be equipped with an electric back angle adjustment option. This function is activated using the control unit in the “additional function” mode.

With the Recaro® LT back rest, the lateral guide in the lumbar section can be individually adjusted by turning the hand wheel. Turning the wheel to the front evenly moves the two lateral guides closer to each other. Turning the wheel to the rear adjusts the two lateral guides away from each other.

Seat

The type X and W seats have an extendable pad at the front. The lever used to adjust this pad is located below the seat. Pull the lever up to shift the pad into the desired position. The pad locks into place as soon as the adjustment lever is released. The adjustment lever must audibly lock in place and return to its initial position.

Headrest

The headrest is adjustable in both height and angle.

A second person is needed to remove the headrest. The pressure points must be detected under the cover by feel and pressed simultaneously (Figure 34). Now the second person can pull the headrest up and out.

Figure 34 Illustration of pressure points

Removing the Recaro® Seat from the Chassis

- Pull ahead on the release strap beneath the front part of the seat.
- Loosen the right and left bolts that lock the seat to the frame.
- To release the posterior seat adapter from the retainer bushings, slightly tilt the seat to the rear and slide it back a little.
- Remove the Recaro® seat.

Attaching the Recaro® Seat to the Chassis

⚠ CAUTION

Risk of injury if the locking bolts are not engaged. When mounting the Recaro® seat, the seat adapter must be securely attached to the retainer bushings. The anterior locking bolts have to snap into place up to the key ring.

- Place the posterior edge of the Recaro® seat on the end of the seat frame.
- Slide the seat forward until the posterior seat adapter hangs in the retainer bushings.
- Tilt the seat ahead until the locking bolts engage with the seat frame up to the key ring.

7.2 Caster Wheel Suspension

The B600 power wheelchair is equipped with drive wheel suspension as standard equipment. A front wheel suspension option is available to increase driving comfort. A spring / damping element on each of the two caster wheels improves driving comfort, especially when driving over uneven surfaces, and improves traction (Figure 35).

Figure 35 Caster wheel suspension

7.3 Caster Wheel Swivel Lock

To navigate ramps or enter lifts, the wheelchair must be driven in completely straight line. An option that locks the caster wheels in the straight ahead position so that turning is impossible is available for the B600 power wheelchair.

Engaging the caster wheel swivel lock:

- Slightly rotate the folding lever on the anterior frame to the side, allowing it to jump out of the mid-position.
- The bolt below the anterior frame extends and engages in the front caster fork as soon as the wheelchair reaches the straight ahead position (see Figure 36, left photo).

Now the power wheelchair can only drive straight forward or backward.

Releasing the caster wheel swivel lock:

- Turn the folding lever back to its mid-position; the front caster fork is released so it can move freely.
- The lever engages in the centre of the anterior frame in an unlocked position (see Figure 36, right photo).

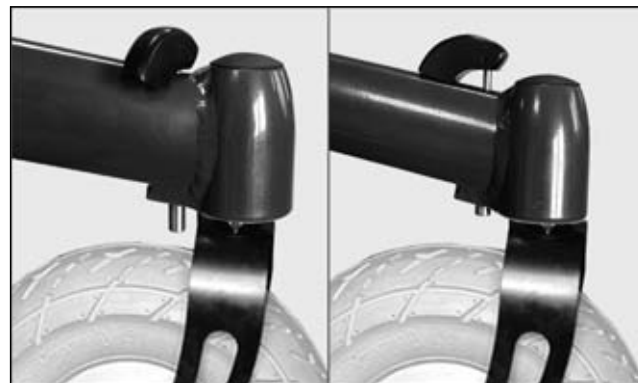
Figure 36 Caster wheel swivel lock engaged and released

7.4 Electrically Adjustable Footrests

INFORMATION

The seat function actuators are not intended for continuous use, but only for short-term limited operation. General recommendation: At maximum load capacity, 10 seconds of operation should be followed by 90 seconds of idle time. The electric seat functions are considered independently of the driving function for this purpose.

To avoid permanent pressure loads or to ensure anti-shock positioning, the power wheelchair offers the possibility of mounting electrically adjustable footrests (Figure 37). With the control unit in “additional function” mode, the right and left footrest can be selected separately or simultaneously. To adjust them, use the joystick.

The amount of space available for getting into or out of the wheelchair can be increased by flipping the footplates up.

Figure 37 Elevated footrest

7.5 Mechanically Elevating Footrest

The B600 power wheelchair can be equipped with mechanically elevating footrests (see Figure 38).

To elevate the footrest:

- Activate the release lever on the footrest (the gas compression spring is activated, see arrow in Figure 38).
- Move the footrest into the desired position.
- Allow the release lever to engage (the gas compression spring is deactivated).

Figure 38 Mechanically elevating footrest

7.6 Control Panel Holder

7.6.1 Height Adjustable Control Panel Holder

An optional special control panel holder can be mounted on the power wheelchair, allowing the height of the control panel to be adjusted. In order to do so, loosen the attachment screw (see Figure 39) and move the control panel to the desired height.

Figure 39 Height adjustment of control panel holder

7.6.2 Swing-Away Control Panel Holder

To allow the user to drive the power wheelchair closer to an object or under the edge of a table, the control panel can be swung to the side with a special control panel holder (see Figure 40).

- Apply slight pressure to move the control panel holder to the side and release the pivot element.
- Swing the control panel holder to the side.
- When bringing the control panel holder back to its original position, the pivot element automatically locks into place.

Figure 40 Swing-away control panel holder

7.7 Special Controls

Various special controls such as the suck-and-blow control, chin control, or foot control can be added to the B600 power wheelchair.

INFORMATION Detailed information is contained in the “Special Controls” instructions for use.

7.8 Additional Indicator and Control Elements

7.8.1 Separate LCD Monitor with Infrared

The separate LCD monitor with infrared allows all types of infrared devices to be controlled using the joystick.

Standard household and HiFi devices can be saved in the control unit using a learning mode. Switches or dimmers are offered as accessories. The mouse emulator is another available infrared module. It allows joystick signals to be used in order to control a computer mouse.

The input devices of the special controls (e. g. sip-and-puff control, chin control) can also be used for infrared functions.

INFORMATION

More detailed information and operating instructions are contained in a separate manual.

Figure 41 LCD module with infrared

Display	Function
	Driving menu with speed level and battery capacity, left direction indicator light activated
	Driving menu with speed level and low battery capacity

Display	Function
	Charging process with drive lock
	Creep speed
	Electric seat height adjustment (optional)
	Electric back adjustment (optional)
	Electric seat tilt (optional)

Display	Function
	Special functions, e. g. coupled electric seat adjustment (back adjustment and seat tilt)
	Electric footrest, left (optional)
	Electric footrest, right (optional)
	Electric footrests, coupled (optional)
	Drive-away lock

Table 6 Information display on the LCD monitor

7.8.2 External Kilometre Indicator

The external kilometre indicator (Figure 42) is attached to the protective bow for the control panel. The kilometre indicator shows:

- Speed
- Trip distance
- Total distance
- Digital clock

INFORMATION When using a separate LCD monitor, an external kilometre indicator is not needed, because the LCD monitor has its own kilometre indicator.

Figure 42 Kilometre indicator

Whenever the display shows km/h or m/h in the upper left corner, the speedometer is activated. Pressing the yellow button lets you see the trip mileage, total mileage, digital clock, and speed one after the other.

Setting the Trip Mileage Back to Zero

Three flashing arrows on the left side of the display indicate the trip mileage counter. To set the value back to zero, press the yellow button for at least 2 seconds.

Setting the Time

If the digital clock is activated on the display, you can set the time. Press the yellow button for at least 2 seconds. The display will flash and alternate between 12: and 24:. Pressing the button again when the display shows 24: will programme the clock to show a 24-hour display. Pressing the button when the display shows 12: will result in a 12-hour display.

The hours and minutes can be set next. The time digits automatically advance one after the other. Press the yellow button when the desired digits are showed on the display.

Setting the Total Mileage Counter

When the total mileage counter is shown on the display, you can switch between km/h and m/h and enter the wheel circumference. To do so, press the yellow button for at least 2

seconds and then press the button again to select between the alternating displays (km/h and m/h). The wheel circumference can be calculated on the basis of the tyre size and can be set in the same way as the clock.

Example for a tyre size of 14 inches:

Conversion from inches to cm	14" x 25.4 = 355.6 mm
C = d x p	355.6 mm x 3.1416 = 1 117 mm (wheel circumference to be set)

7.8.3 Push Button Module

CAUTION

Risk of accidents due to uncontrolled driving characteristics. The push button module is an option for the use of the enAble50 system. Improper adaptation of the push button module to a power wheelchair with components outside the Otto Bock modular system or subsequent modifications by the user are not permitted.

The push button module lets the user select additional electrical functions directly during normal driving operation.

Available functions

Depending on the power wheelchair model and the programming of the push button module, up to 5 electric functions can be controlled during normal driving operation (Figure 43, item 1). Examples include:

- Seat tilt
- Electric back angle adjustment
- Lift seat function
- Right / left footrest; combination of both footrests
- Special functions (combinations) => symbols S1 – S5

The symbols used on the push button module correspond to the symbols used on the LCD monitor (see 7.8.1 "Separate LCD Monitor with Infrared").

Push button functions

The "M" button (= mode, Figure 43, item 2) is used to advance the individual functions (1-2-3-4-5-1-2-...).

The blue LED (Figure 43, item 3) indicates the selected function.

The selected function is executed with the arrow keys (Figure 43, item 4) (function up / down).

INFORMATION

The function of the arrow keys can be adapted according to the wishes of the user (up arrow = function up or function down – depending on programming). Programming may only be completed by instructed specialist staff.

Three freely positionable Buddy buttons may also be connected to the push button module as accessories (Figure 43, item 5). They assume the function of the mode and up / down arrow keys (Figure 43, items 2/4). Symbols on the back of the push button module indicate the corresponding function.

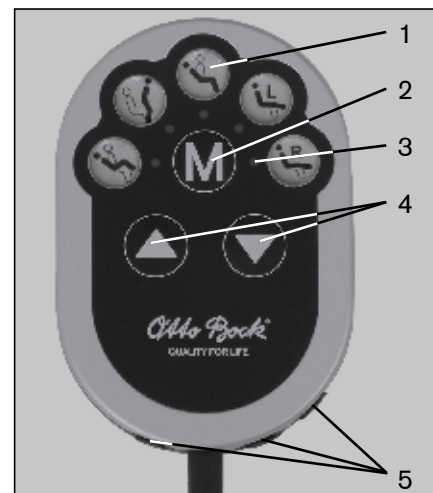
Figure 43 Push button module

- 1 Available electric functions
- 2 "Function selection" mode button
- 3 "Selected function" LED indicator
- 4 Up / Down arrow keys (see information box)
- 5 Stereo plug for the connection of Buddy buttons

7.8.4 2x Adapter Cable for Piko or Buddy Button

The power wheelchair may be equipped with a 2x adapter cable. The respective Piko or Buddy button that is connected may be used to control electrical seat functions. The functions to be controlled are selected when the power wheelchair is ordered. Your specialist dealer can change the assignment of the button using a handheld programming device or programming station.

The buttons are equipped with a so-called toggle function. This means that the selected function is activated as long as the button is pressed.

Pressing the button again activates the selected function in the opposite direction. For example, when the function for the electric back rest is activated, the back rest moves back as long as the button is pushed the first time the button is activated. Activating this function again causes the back rest to move forward as long as the button is pushed.

7.9 Additional Options

7.9.1 Lap Belt

WARNING

Risk of injury due to improper use of the lap belt. The lap belt helps to additionally stabilise the person sitting in the power wheelchair. The lap belt must never be used as part of a restraint system for transportation of the person in a motor vehicle for the disabled.

The B600 can be equipped with a lap belt. The lap belt serves exclusively as a safety belt when driving the power wheelchair.

To secure the lap belt, push the two halves of the buckle together so they engage (Figure 44). The belt buckle must snap in place audibly. Verify that the belt has locked by trying to pulling it apart. The lap belt should not be too tight on the body. Objects caught under the belt can cause painful pressure sores.

Release the lap belt by pressing the red unlock button.

The belt length can be adjusted on both sides. To set the length, position both buckle halves in the middle of the body. After positioning the buckle halves at a right angle in relation

to the belt, their position can be varied. Any excessive belt length will be held in place by the plastic slides.

Figure 44 Applying the lap belt

7.9.2 Four-Way Chest Strap

⚠ WARNING

Risk of injury due to improper use of the four-way chest strap. The four-way chest strap helps to additionally stabilise the person sitting in the power wheelchair. The four-way chest strap must never be used as part of a restraint system for transportation of the person in a motor vehicle for the disabled.

The four-way chest strap makes it possible to restrain the patient in the power wheelchair. The lower strap should be placed between the pelvic bone and thigh, and should not be too tight. The belt buckle is positioned at the front of the body, in the centre. The shoulder straps should be positioned evenly on both sides on the thorax, and should not be too tight.

To close the four-way chest strap, connect all parts of the belt buckle. The belt buckle must snap in place audibly. Verify that the belt has locked by trying to pulling it apart. The four-way chest strap should not be too tight on the body. Objects caught under the belt can cause painful pressure sores.

To increase the length of the belt, turn the buckle by 90° and

pull. To reduce the length of the belt, pull the overhanging belt end.

7.9.3 Adapter for Headrest Mounting Kit (Standard Seat)

The power wheelchair can be equipped with an optional headrest. An adapter for attachment of the headrest mounting kit is available for this purpose.

Figure 45 Adapter for headrest mounting kit

7.9.4 Attendant Control

A separate control panel can be attached to the back rest of the power wheelchair for use by an attendant.

7.9.5 Additional Options

- Curb climbing assist: The curb climbing assist is used to clear curbs and steps that have a maximum height of 10 cm.
- Puncture resistant tyres: Solid rubber tyres
- Contour seats in various sizes for improved comfort and additional adjustment options
- Seating shell interface: For the adaptation of special seating shells. Control panel holders for seating shells are also available.
- Arm rest accessories: Special adapter for arm rests from our accessory catalogue
- Attendant control: External joystick on the back rest
- Joystick accessories: Fork for tetraplegics, STICK S80, soft ball, golf ball, and flexible control stick shaft
- Protective bow for control panel: Metal bow to protect the control panel against impacts
- Crutch holder
- Folding rear view mirror
- Tray: Attachable table tray
- Accessories case, backpack, pocket for mobile phone

These and other optional add-on components are included on the order form and in the wheelchair accessories catalogue.

8 Malfunctions / Troubleshooting

INFORMATION

Should you encounter problems while troubleshooting or if you do not manage to completely eliminate a problem by following the measures described here, please contact your specialist dealer.

Malfunctions are indicated on the control panel LCD monitor. Table 7 describes the individual error codes with the corresponding problem sources as well as the possible causes and resolutions.

If malfunctions cannot be fully corrected using the procedures described here, an authorized dealer has the ability to read the exact error code with a hand-held programming device and to carry out a targeted system analysis.

All problems that have ever occurred are recorded and can be retrieved, e.g. in case of a general overhaul of the power wheelchair. For example, the saved data can be used to determine future service and maintenance intervals.

8.1 Warning

A warning indicates a status or malfunction of one or several components of the power wheelchair. The function of components that have no errors is not restricted. For example, if the connection between the controller and seat motor malfunctions, this error will only be indicated if the user attempts to activate the motor. However, the driving function is still available.

8.2 Error

⚠ CAUTION

Risk of injury due to stopping the wheelchair abruptly.


In the event of communication problems in the bus system of the control unit, the system triggers an emergency stop to prevent any uncontrolled functions. Depending on the source of the malfunction, it is possible to drive the wheelchair out of a danger zone such as road traffic after switching the control unit on again.

If the driving function is still not available after switching the control unit on again, the brake must be deactivated (see Section 6.4) to switch over to pushing mode. After the power wheelchair has triggered an emergency stop, contact a specialist dealer as soon as possible!

⚠ CAUTION

Risk of accidents due to uncontrolled driving characteristics. Uncontrolled movements can occur during the operation of the power wheelchair as a result of malfunctions. In this case, please contact your authorised specialist dealer immediately.

An error affects one or several functions of the power wheelchair. The system is not fully operational until the error has been corrected.

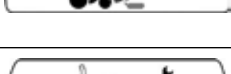Display symbol(s)	Error / warning	Cause	Possible Corrective Action
	Controller temperature warning	Overheating due to excessive load	Cool down phase
	Motor temperature warning	Overheating due to excessive load	Cool down phase
	Joystick warning	Joystick not in zero position when switching on	Bring joystick to zero position before switching on
	Manual control device fault	Defective joystick	Contact specialist dealer
	Controller error	Defective controller	Contact specialist dealer
	Back angle adjustment motor fault	Faulty cabling or plug connection; defective actuator	Check cabling / plug connection; contact specialist dealer
	Lift seat motor fault	Faulty cabling or plug connection; defective actuator	Check cabling / plug connection; contact specialist dealer
	Seat tilt motor fault	Faulty seat tilt function, cabling, or plug connection; defective actuator	Check cabling / plug connection; contact specialist dealer

Display symbol(s)	Error / warning	Cause	Possible Corrective Action
	Electric footrest motor fault	Faulty cabling or plug connection; defective actuator	Check cabling / plug connection; contact specialist dealer
	Drive motor fault	Faulty cabling or plug connection; defective actuator	Check cabling / plug connection; contact specialist dealer
	Brake fault	Open brake release / defective brake	Close brake release; check brake (e.g. the cable)
	Emergency stop	Severe fault caused by controller, manual control device, and / or drive motor malfunction	Contact specialist dealer
	Communication error (alternating flashing signal)	Faulty connection between manual control device and controller / defective cabling, software or hardware	Check cabling / plug connection; contact specialist dealer
	Battery under voltage	Battery deep discharge	Charge as soon as possible
	Battery over voltage	Voltage too high (after full charge and driving downhill)	Continue driving slowly

Table 7 Status and error messages

8.3 Defect / Failure

A defect represents a serious failure of a system component. A failure is the most serious type of malfunction, which causes an immediate emergency stop of the system.

Defects / failures are shown by a continuous sequential light indicator with an audible signal.

After the malfunction has been corrected, the system is re-activated by restarting it.

9 Maintenance and Care

INFORMATION

A spare parts catalogue from Otto Bock is available upon request for ordering spare parts. Only spare parts supplied by Otto Bock may be used. Failure to comply with this requirement will render the warranty null and void.

INFORMATION

Should you encounter problems during maintenance, contact your specialist dealer. The power wheelchair is to be checked and serviced once a year by a specialist dealer.

9.1 Maintenance Intervals

The correct functionality of the power wheelchair must be checked prior to each use. The items listed in Table 8 must be checked by the user at the indicated intervals.

Component	Activity	Daily	Weekly	Monthly
Arm rest and side panel	Fastening screws tightened Arm rest and control panel secured Check arm rest for damage		Prior to every use X	X
Drive wheels	Wheels must rotate freely and without axial run out Central nut on the driving shaft tightened Check if wheel mounts are seated securely Directional stability of the entire wheelchair		X	X X X
Tyres	Air pressure (printed on the sidewall of the tyre) Sufficient tread depth, at least 1 mm Check for damage			X X X
Batteries	Check fluid / acid level (not with gel batteries)			X
Lighting	Visually check for damage Test for functionality	X	X	
Electronics	Control unit / push button module free of error messages Charger does not indicate any error messages on the LCD monitor Check plug connections		Prior to every use X	X

Component	Activity	Daily	Weekly	Monthly
Brake	Activate brake lever while control unit is switched on Brake function active with brake engaged	X		X
Mechanically elevating footrest	Check ratchet mechanism for functionality and secure seating Check footplates for damage Visual check for scratches on the piston rod and oil leakage			X X X
Lift seat	Visual inspection of all moving components, especially cabling - check for damage Check whether screw connections are tight			X X
Caster wheels	Fork seated in the adapter without play Wheels must rotate freely and without axial run out Fastening nut tightened			X X X
Padding and belts	Proper condition of padding No wear on the seat belts Check belt buckle for functionality		X	X X
Seat attachment	Check if attachment screws are seated securely Check seat lock of Recaro® seat			X X

Table 8 Maintenance measures and intervals

9.2 Changing the Fuse

The 80 A fuse is located in a fuse holder at the rear end of the battery case (see Figure 46).

Open the cap of the safety holder and exchange the fuse. Make sure that the fuse is pressed into the centre of the spring contacts and that it is not positioned at an angle. Close the cap until you feel it snap back into place.

Figure 46 Fuse holder

- 1 Cover open
- 2 Fuse inserted
- 3 Fuse holder

9.3 Changing Tyres

INFORMATION

Direct sunlight (UV light) causes the tyres to age prematurely. As a result, the tread surface hardens and corner pieces break out of the tread.

INFORMATION

Avoid parking the wheelchair outdoors whenever possible. The tyres should be replaced at least every 2 years regardless of wear and tear.

If the wheelchair is parked for an extended period of time or the tyres overheat (e.g. in the vicinity of radiators or in case of exposure to strong sunlight behind glass), the tyres may become permanently deformed. Therefore, sufficient clearance should always be maintained from sources of heat, and the wheelchair should be moved frequently or placed on blocks for storage.

To change the tyre of a drive wheel, proceed as follows:

- Secure the wheelchair against tipping sideways by using a suitable block under the drive unit sustainer.

- To remove a drive wheel, loosen the 4 Allen head screws in the middle of the wheel (Figure 47) and remove the wheel.
- If the inner tube of a drive wheel needs to be replaced, loosen all Allen head screws on the inside of the rim and pull apart the two-piece rim.

The defective inner tube is now freely accessible and can be replaced.

Figure 47 Disassembling the drive wheel

To change the tyre of a caster wheel, proceed as follows:

- To remove the caster wheel, loosen the axle screw with a size 6 Allen key (see Figure 48, item 1) and pull out the axle.
- Loosen all Allen head screws (see Figure 48, item 2) and separate the two-piece rim.

The defective inner tube is now freely accessible and can be replaced.

Figure 48 Disassembling the caster wheel

- 1 Axle screw
- 2 Allen head screw

9.4 Changing Defective Lights

NOTICE

Risk of damage due to moisture. In order to prevent moisture from penetrating the lamp, ensure that the pane is seated properly in the housing during installation and that the screws are tightened.

INFORMATION

Lamp housings or lamps can be ordered from a specialist dealer.

To change the halogen lamp in the front light, proceed as follows:

- Fold the pane forward with slight pressure on the black locking lever at the lower end of the lamp (see Figure 49, item 1).
- Replace the defective halogen lamp with a new one.
- Close the pane.

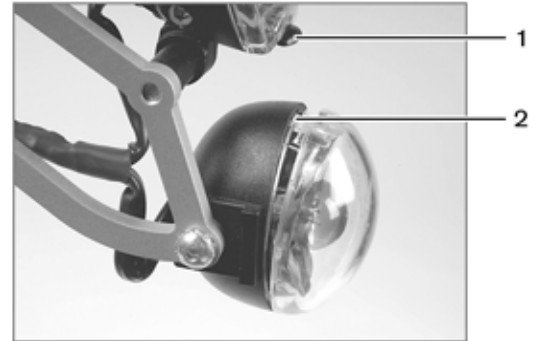
Figure 49 Changing the front lights

- 1 Front light locking lever
- 2 Recess (use the screwdriver at this point)

To change the lamp of the front direction indicator, proceed as follows:

- Open the front direction indicator: Insert a narrow screwdriver into the recess on the direction indicator housing and tip the pane down (see Figure 49, item 2).
- Remove the pane.
- Slightly rotate the lamp to release it from the lock on the lamp socket and pull it out (see Figure 50).

- To install the lamp, slide it in and rotate it into the lock.
- Insert the lamp socket into the housing and lock the pane of the direction indicator in place

Figure 50 Replacing the direction indicator lamp

To change the rear direction indicators and the rear light, proceed as follows (see Figure 51):

- Unscrew the pane.
- Press the lamp in slightly and give it a quarter turn counter-clockwise, then pull it out.
- Insert new lamp and give it a quarter turn clockwise.

Figure 51 Rear lamp, pane removed

9.5 Cleaning and Care

NOTICE

Risk of damage to electronics due to water penetration. In order to avoid malfunctions, the electronic components, motor, and batteries must not come into contact with water when cleaning the power wheelchair.

NOTICE

Risk of damage to power wheelchair components. Do not use any aggressive cleaners, solvents, or hard brushes to clean the power wheelchair. Never use a water jet or high-pressure cleaning apparatus to clean the power wheelchair.

INFORMATION

Prior to disinfection, clean the seat and back upholstery as well as the seat cushion, control panel, and arm rest.

The power wheelchair must be cleaned regularly, depending on the amount of use and the degree of soiling.

The control panel, battery charger, arm rests, and trim can be cleaned with a damp cloth and mild cleaning solution.

Use a dry brush to clean the seat and back upholstery as well as the seat cushion.

Use a damp plastic brush to clean the wheels and frame.

10 Technical Data

Dimensions and Weights	
Seat width	38 – 42 cm or 43 – 48 cm
Seat depth	38 – 46 cm or 42 – 50 cm
Seat height	45 – 60 cm
Arm rest height	24 – 36 cm
Arm rest length	26 cm
Lower leg length	25 – 34 cm or 35 – 44 cm or 45 – 54 cm
Back height	45 oder 55 cm
Back angle	-9/1/11/21° or 0/10/20/30°
Overall width	64.5 cm
Overall height	103 cm
Total length	108 cm
Turning radius	155 cm
Tyre size	
Steering caster	9"-10"
Drive wheel	14"
Air pressure	Front: Printed on the sidewall of the tyre
Weight empty	115 kg (depending on options)
Maximum load capacity	140 kg (optionally 180 kg)

Electrical Equipment	
Operating voltage	24 V
Batteries:	
Wet cell batteries	2 x 12 V, 60 Ah
Gel batteries	2 x 12 V, 63 Ah (C5)
Control unit:	
Model	enAble50
Operating voltage	24 V DC
Max. output current per motor	130 A
Push button module:	
Model	For enAble50 control unit
Protection class	IP 64
Operating and storage temperature range	-20°C to + 80°C
Lighting:	
Front direction indicator	H21W 12 V BAY9s
Front light	HMP 08 2,4 W; 6 V, PX13,5s
Rear direction indicator	C21W 12 V; BA15s
Rear light	C5W 6 V; BA15s
Fuse	80 A

Driving Data	
Speed	6 km/h or 10 km/h
Climbing ability	17 %
Maximum obstacle height	5 cm
Range	Approx. 35 km
Operating temperature	-25 °C to +50 °C
Transport and storage temperature range	-40 °C to +65 °C
Battery Charger	
Model	MEG 2409 automatic battery charger with computer-controlled characteristic curve
Power requirements	230 V -10 %, +6 %
Mains frequency	50 Hz ±4 %
Protection class	1 (protective conductor)
Charging connection	24 V DC
Nominal charging current	9 A
Residual ripple	< 1%
Characteristic curve	IUU with trickle charging phase, computer-controlled
Primary fuse	T 4 A fuse link, not accessible from outside

Battery Charger	
Secondary circuit breaker	Electronic, reversible reverse battery protection, short circuit-proof, idling-proof, overheating protection
Protection class	IP 21
Ambient temperature	-10 °C to +40 °C
Display	2 LED
Weight	1.7 kg
Dimensions (WxHxD)	105 x 65 x 205 mm

Table 9 Technical data

11 Disposal

INFORMATION

If a wheelchair is to be disposed of, all components and materials of the power wheelchair must be recycled or disposed of properly.

⚠ CAUTION

Risk of polluting the environment with battery acid.

The batteries of the power wheelchair contain a toxic acid. They must not be disposed of with regular domestic waste and the battery acid must not enter the sewage water system or ground. You must observe the battery manufacturer's instructions printed on the batteries.

If the power wheelchair is no longer in use, it must be disposed of properly in accordance with national environmental protection regulations.

Please return defective batteries to your specialist dealer when buying new ones.

12 Information on Re-Use

The B600 power wheelchair is suitable for re-use.

Similar to second-hand machines or vehicles, products that are being re-used are subject to increased strain. Features and functions must not change in a way that could endanger patients or other persons within the product's life cycle.

Based on market observations and the current state-of-the-art, the manufacturer has calculated that the B600 power wheelchair can be used for a period of **5 years**, provided that it is used properly and that the service and maintenance instructions are observed. Periods during which the wheelchair is stored at the dealer or with paying parties are not included in this period. It should be clearly pointed out, however, that the B600 is a reliable product far beyond this defined period of time, provided that it is cared for and maintained properly.

In cases of re-use, the corresponding product must first be thoroughly cleaned and disinfected.

Subsequently, the product must be examined by an authorized specialist to check the condition, wear, and possible defects.

Any worn and damaged components as well as components which do not fit or are unsuitable for the new user must be replaced.

The service manual includes a service schedule for each model, detailed information, and a list of the required tools.

13 Liability

The manufacturer's warranty applies only if the device has been used under the conditions and for the purposes described. The manufacturer recommends that the device be used and maintained according to the instructions for use.

The manufacturer is not responsible for damages caused by components and spare parts not approved by the manufacturer. Repairs must be carried out exclusively by authorised dealers or by the manufacturer.

14 CE Conformity

This device meets the requirements of the 93/42/EEC guidelines for medical devices. This device has been classified as a class I device according to the classification criteria outlined in appendix IX of the guidelines. The declaration of conformity was therefore created by Otto Bock with sole responsibility according to appendix VII of the guidelines.

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**Otto Bock HealthCare
Canada Ltd.**
CA-Canada, Burlington,
Ontario L7L 5N5
Tel. +1 289 288-4848
Fax +1 289 288-4837
info@ottobock.ca

**Otto Bock HealthCare
Andina Ltda.**
Bogotá / Colombia
Tel. +57 1 8619988
Fax +57 1 8619977
ottobock@telesat.com.co

**Otto Bock de Mexico
S.A. de C.V.**
MEX-Guadalajara, Jal. 44210
Tel. +52 33 38246787
Fax +52 33 38531935
info@ottobock.com.mx

Otto Bock HealthCare LP
U.S.A.-Minneapolis,
Minnesota 55447
Phone +1 800 328 4058
Fax +1 800 962 2549
usa.customerservice@
ottobock.com

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Otto Bock Mobility Solutions GmbH has been certified by the German Society for the Certification of Quality Assurance Systems (DQS) in accordance with DIN EN ISO 13 485 : 2007, reg. no. 435514 MP27 (management system)

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HealthCare Inc.**
South East Asia Co. Ltd.**
Otto Bock HealthCare GmbH
Otto Bock Manufacturing Königsee GmbH